




# RB780 and RY775 can be used together in flow cytometry panels


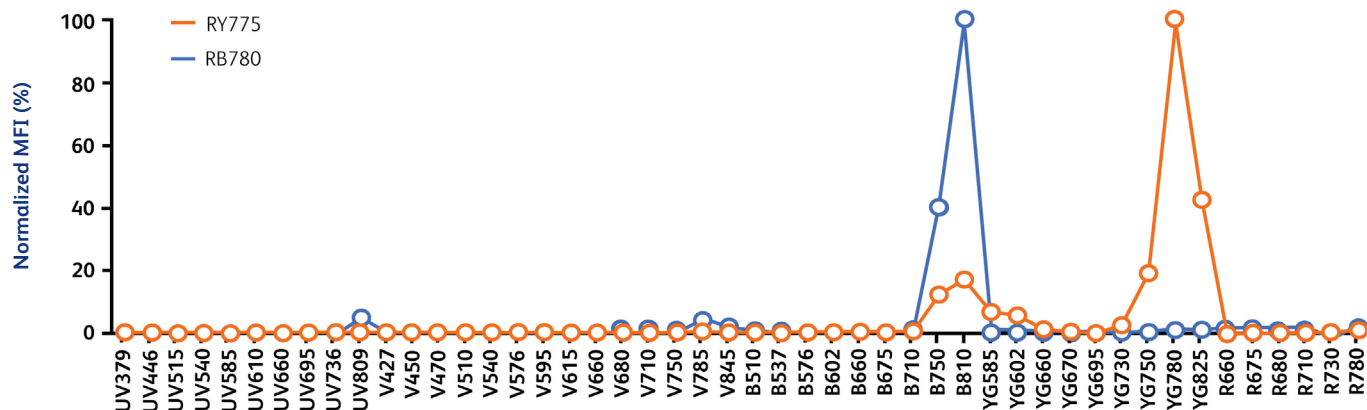


## Normalized emission profiles for RB780 and RY775

RB780 and RY775 have distinct emission profiles and minimal cross-laser excitation. Samples run on a BD FACSymphony™ A5 SE Cell Analyzer.

For improved performance and increased panel design flexibility, RY775 reagents can be paired with BD Horizon RealBlue™ 780 (RB780) Reagents. Both dyes have reduced cross-laser excitation, which enables their superior performance when used together on instruments with both blue and yellow-green lasers and the appropriate filters, such as the BD FACSymphony™ A3, A5 and A5 SE Analyzers as well as the BD FACSDiscover™ S8 Cell Sorter.

## Pair RB780 with RY775 instead of PE-Cy7 for more flexibility in panel design

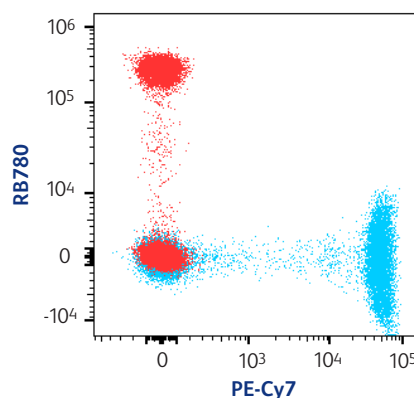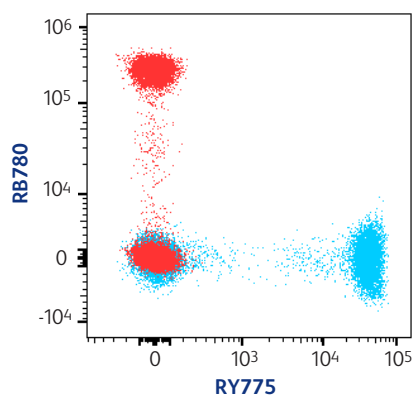

**RY775 improves performance when used with RB780 on instruments equipped with both blue and yellow-green lasers**

PBMCs stained with either CD4 RB780, CD4 RY775 or CD4 PE-Cy7. Samples acquired on the BD FACSDiscover™ S8 Cell Sorter and spectrally unmixed in FlowJo™ Software with BD SpectralFX™ Technology.



For more information, please visit [bdbiosciences.com/real](https://bdbiosciences.com/real).

BD flow cytometers are Class 1 Laser Products.

For Research Use Only. Not for use in diagnostic or therapeutic procedures.

BD, the BD Logo, BD FACSDiscover, BD FACSymphony, BD Horizon RealBlue, BD Horizon RealYellow, FACS, FlowJo, Pharmingen and SpectralFX are trademarks of Becton, Dickinson and Company or its affiliates. All other trademarks are the property of their respective owners. © 2024 BD. All rights reserved. BD-121944 (v1.0) 0424

Cy is a trademark of Global Life Sciences Solutions Germany GmbH or an affiliate doing business as Cytiva.

